



# Finally... A Safe, Simple Way to Remove Scalpel Blades

**OSHA**  
*Compliant*



 **QLICKSMART**<sup>®</sup>  
BLADE REMOVAL SYSTEMS

**Mopec** The Better Product People

Call Us at: **1-800-362-8491** or Visit Us on the Web at **www.mopec.com**

# The World's First Single-

# For Safe and




- Puncture Resistant
- Spill Proof
- Self Contained Unit
- Maintenance Free
- Built-In Blade Counter Prevents Overfilling
- Accommodates Standard Blades and Blade Handles

*The revised OSHA Bloodborne Pathogens Standard 29 CFR Part 1910.1030 mandates that:*

- (d)(2)(vii) Contaminated needles and other contaminated sharps shall not be bent, recapped, or removed except as noted in paragraphs (d)(2)(vii)(A) and (d)(2)(vii)(B) below. Shearing or breaking of contaminated needles is prohibited.
- (d)(2)(vii)(A) Contaminated needles and other contaminated sharps shall not be bent, recapped or removed unless the employer can demonstrate that no alternative is feasible or that such option is required by a specific medical or dental procedure.
- (d)(2)(vii)(B) Such bending, recapping or needle removal must be accomplished through the use of a mechanical device or a one-handed technique.
- (d)(2)(viii) Immediately or as soon as possible after use, contaminated reusable sharps shall be placed in appropriate containers until properly reprocessed. These containers shall be:
  - (d)(2)(viii)(A) Puncture resistant;
  - (d)(2)(viii)(B) Labeled or color-coded in accordance with this standard;
  - (d)(2)(viii)(C) Leakproof on the sides and bottom; and
  - (d)(2)(viii)(D) In accordance with the requirements set forth in paragraph (d)(4)(ii)(E) for reusable sharps.
- (d)(4)(iii)(A)(1) Contaminated sharps shall be discarded immediately or as soon as feasible in containers.
- (d)(4)(iii)(A)(2)(i) Easily accessible to personnel and located as close as is feasible to the immediate area where sharps are used or can be reasonably anticipated to be found (e.g., laundries);
- (d)(4)(iii)(A)(2)(iii) replaced routinely and not be allowed to overfill.

# Handed Scalpel Blade Remover

## Easy Blade Disposal

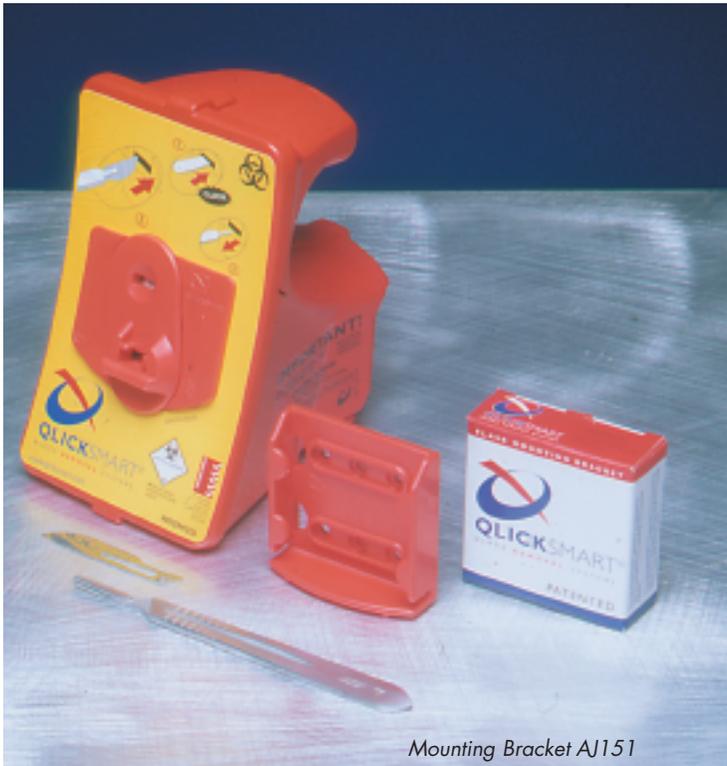

### Finally... A Safe and Simple Way to Remove Scalpel Blades

The Clicksmart Blade Removal System is a medical safety device specially designed, engineered and manufactured to comply with recently revised OSHA Bloodborne Pathogens Standards.

Clicksmart offers single-handed removal and immediate containment of contaminated scalpel blades which dramatically reduces the risk and cost of injury to the user.

### Simply Amazing

Clicksmart is safe, inexpensive and easy to use. Simply insert the scalpel into the slot, listen for the "click" and then withdraw the handle. The blade is securely contained inside the sealed, spill proof and puncture resistant container.

### Safe and Easy Disposal

An automatic shut-off feature prevents the Clicksmart system from overfilling and a built-in counter clearly displays the remaining capacity (100 blades). When full, simply dispose of the fully sealed self contained unit as medical waste.

### No Special Blades or Handles

Clicksmart has been designed to work with most standard scalpel instruments, whether they are Mopec's or a competitor's, including blades 10, 10A, 11, 12, 15, 20, 21, 22, 23, 24 and handles 3, 4, 4L, 7, 9.

### Counter, Cart or Wall Mount

The universal mounting bracket and positioning template allows quick and convenient counter, cart or wall-mount. The industrial-strength mounting bracket adhesive provides immediate and simple installation. Supplied with optional screws and plugs.

### Maintenance Free

Fully self-contained and sealed Clicksmart has no mechanical device to maintain, break or fix. Simply replace the unit when full.

**Clicksmart Product No. AJ150**

**Mounting Bracket Product No. AJ151**


**Fully self-contained and sealed, Qlicksmart has no mechanical device to maintain, break or fix. Simply replace the unit when full.**



**The World's First Single-Handed Scalpel Blade Remover with the added benefits:**



*Check out our Catalog for a Complete Selection of Blades and Handles That Will Fit the QlickSmart Blade Removal System*

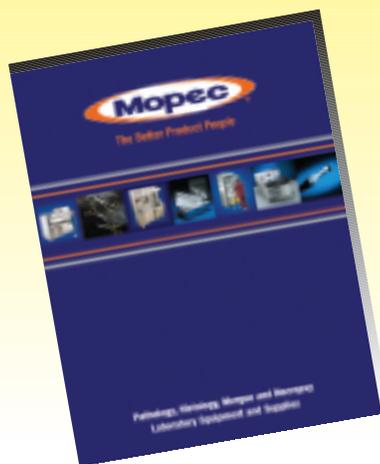

## Now Available - The New Mopec Catalog

### Inside you'll find:

- Over a hundred new products
- Complete descriptions of all our products
- New, ergonomic workstations and tables
- Products and workstations you'll only find at Mopec

For a complete catalog detailing our products and services, please call or visit our web site



21750 Coolidge Highway • Oak Park, MI 48237  
Toll Free: 800-362-8491 • Tel: 248-291-2040 • Fax: 248-291-2050

[www.mopec.com](http://www.mopec.com)